




XYLOTRIOSE (Lot 180803)

O-XTR

10/18

MW: 414.36

CAS: 47592-59-6

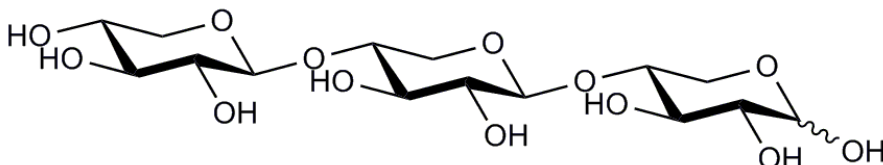

PREPARATION:

Prepared by controlled acid hydrolysis of arabinoxyylan.

PURITY: > 95% (HPLC)

HPLC:


Column: 2 x Tosoh TSK-GEL G2500 PWXL (7.8 x 300 mm) plus guard column (7.8 x 35mm)

Temperature: 80°C

Mobile phase: dH₂O

Flow rate: 0.5 mL/min

HPLC System: Waters Alliance e2695 Separations Module, Waters 2414 RI detector and Empower v 3 software

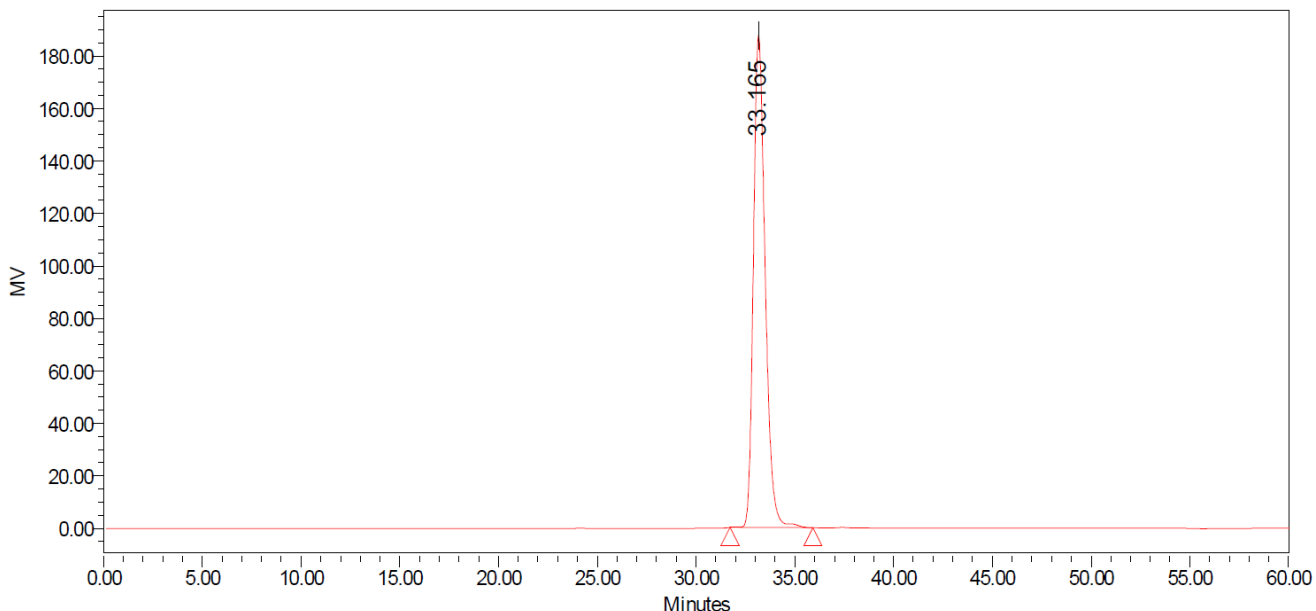

Peak Name:

	Injection	RT	Area	% Area	Height
1	1	33.165	7735867	100.00	187467

HPAEC-PAD:

Column: CarboPac PA200 guard and analytical columns (3 x 250 mm)

Temperature: 30°C

Detector: Au electrode; waveform Carbohydrate, standard quad

Flow rate: 0.5 mL/min

IC system: Dionex ICS5000+ DP system and Chromeleon 7 software

A stepwise linear gradient method was employed as shown.

Time (min)	100 mM NaOH (%)	120 mM NaOAc (%) in 100mM NaOH
0	100	0
10	65	35
14	50	50
15	0	100

